


The Panacea® Pro Cushion

**DIRECT
SUPPLY®**
Equipment & Furnishings

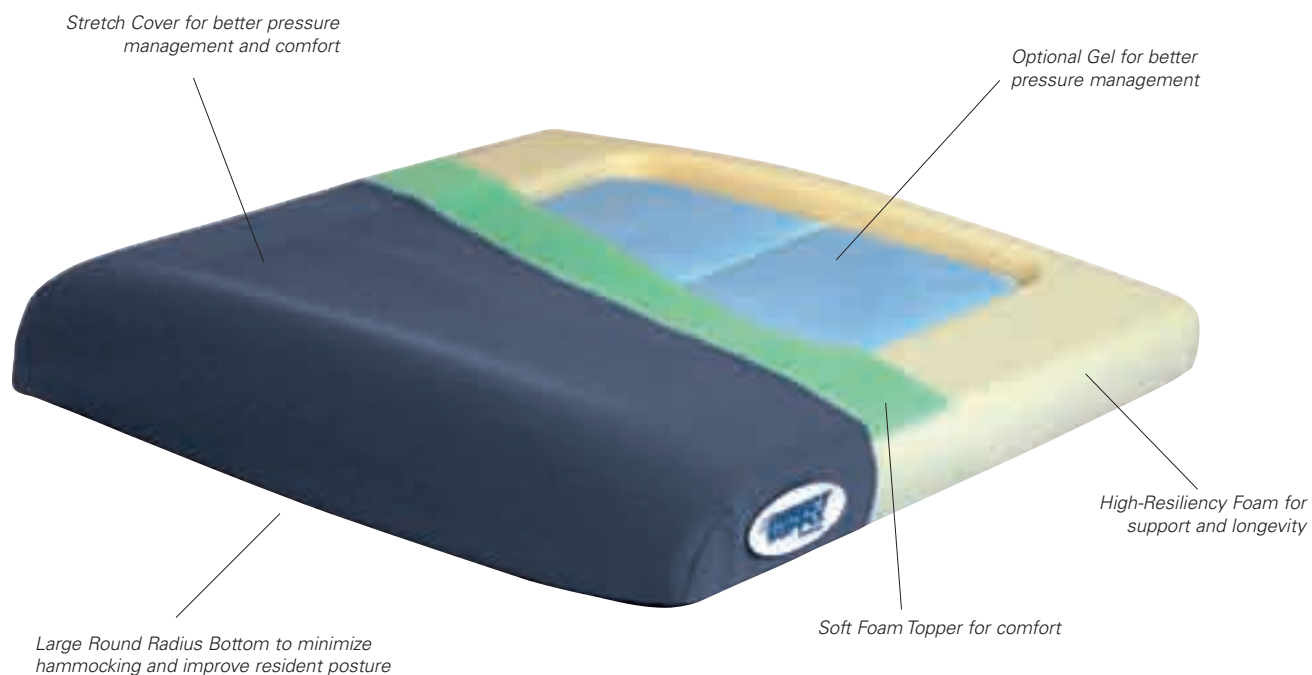

Panacea® Wheelchair Cushions & Positioning Devices

Comfort and Quality at a Great Value


Panacea® Pro Cushions

**DIRECT
SUPPLY®**
Panacea.

Panacea® Air Cushions



Pommel

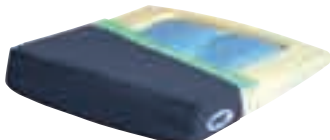

Wedge



Anti-Thrust



Saddle



Adjustable Air



Convertible

1-800-634-7328

directsupply.com

CUSHIONS / POSITIONING

Benefits

Stretch Cover: The stretch cover on Panacea Cushions offers great immersion in the cushion for better pressure redistribution. It also keeps the seating surface cooler and more comfortable, and is moisture resistant. Most cushions include a strap to keep the cushion securely in place.

Optional Gel: Many Panacea Cushions are available with gel to increase surface area contact with the user which improves pressure redistribution and draws heat away from the resident for a cooler seating surface. Optional chambered gel provides added stability and keeps the gel in position where it's most needed. Additionally, gel packs on all Panacea Cushions are placed at the top of the cushion to provide the most benefit, rather than buried beneath thick foam like on other cushions.

Large Round Radius Bottom: Panacea Cushions feature a standard round bottom to fit the sling seat found on most wheelchairs. This minimizes hammocking and keeps the hips better aligned while also eliminating the need for a solid seat insert.

High-Resiliency Foam: Panacea Cushions feature high-quality foam not found on other cushions. It will maintain its integrity and shape, and hold up over years of use. Additionally, it provides the support needed to achieve proper positioning.

Soft Foam Topper: Many Panacea Cushions feature a soft foam topper over the high-resiliency base and optional gel. This topper conforms to the resident's individual curves and provides a more comfortable seating environment.

Molded Foam: Many Panacea Cushions use new molded foam technology to provide improved resistance levels as the user immerses into the foam. This technology has significantly improved cushion comfort and long term durability.

Styles

Flat / Zero Elevation: These cushions are the same thickness across the entire surface. They are designed for residents who need added comfort, support, or basic pressure management with minimal positioning.

Wedge: These cushions are thicker in the front than in the back. They are designed to provide support, comfort, and pressure reduction for active residents who tend to thrust or slide forward. They are an excellent match with reclining wheelchairs to distribute pressure away from the tender ischium and coccyx areas to the more supportive back.

Anti-Thrust: Anti-thrust cushions have a 1" shelf at midline, and are designed to reduce sacral sitting or hip rotation and prevent forward sliding.

Pommel: Pommel cushions abduct and properly align legs which promotes better posture. A common term for incorrectly positioned legs is "wind-swept legs". This can cause the hips to rotate, the spine to twist and the head to drop, creating many long term health problems.

Saddle: Saddle cushions have a geometric design with contour leg abductors and pommel. They provide proper alignment, reduce hip rotation and pelvis obliquity, and increase seating area to reduce pressure.

Panacea® Wheelchair Cushions

Panacea Comfort Memory Foam – Conforming Comfort

Stretch Cover for better pressure management and comfort

Visco-Elastic Memory Foam conforms to residents' contours




High-Resiliency Foam Base for support and longevity

Panacea Bariatric – Comfortable support for residents up to 650 lbs.

Stretch Cover for better pressure management and comfort

Soft, conforming foam or optional gel

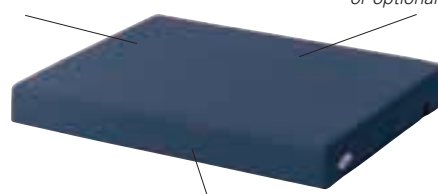


High-Density Foam Base for support and longevity

Panacea® Wheelchair Cushions


Panacea Position Flo – Contoured positioning and pressure management



Panacea Performance – Advanced positioning and pressure management

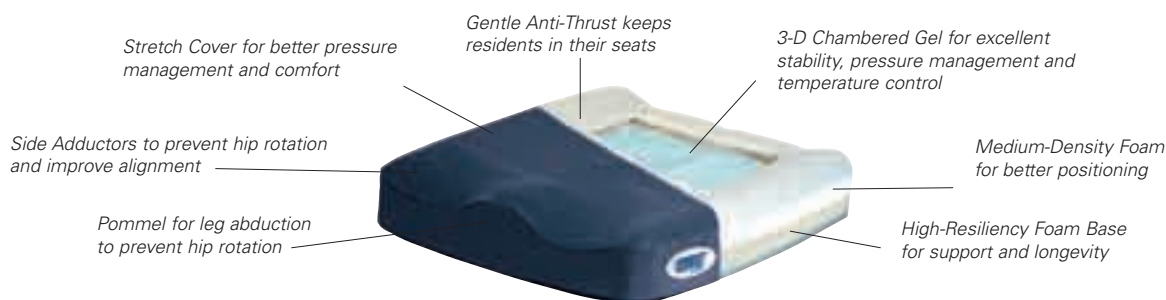


Panacea 4-Layer with Gel

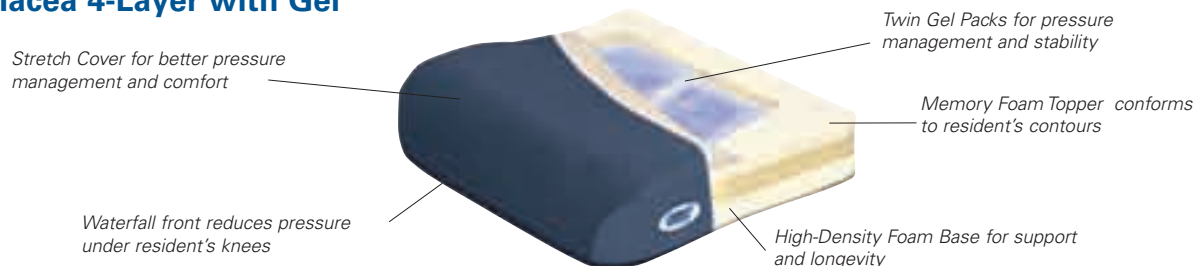

Panacea Adjustable Air Cushion – for advanced pressure management



Panacea Convertible Cushion – with optional alternating pressure

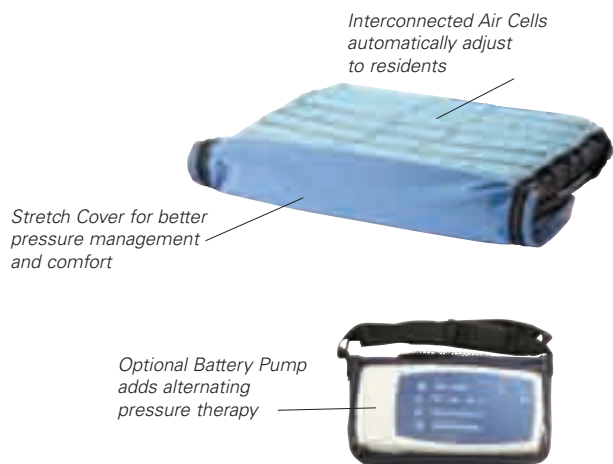



Panacea® Lateral Supports – provide mild to moderate trunk support

- Easy-to-install lateral support
- Pads adjust laterally and vertically with Velcro®
- Available with 6"W or 4"W laterals

#97702 Panacea Lateral Support System



Panacea® Lap Supports – fit all Invacare & Easy Care wheelchairs!

- Soft top cushion fits between resident and chair for upper body positioning
- Provides comfortable elbow rest and helps prevent forward leaning

#84427 Lap Support with Full Length Armrest (for 16"-18"W Wheelchairs)

#84428 Lap Support with Desk Length Armrest (for 16"-18"W Wheelchairs)

#19425 Lap Support with Full Length Armrest (for 22"W Wheelchairs)

#19426 Lap Support with Full Length Armrest (for 24"W Wheelchairs)



Panacea® Arms – anatomically contoured to cradle resident's arm

- Crowned palm area to avoid finger contractures
- Arm rotates 270°
- Soft molded urethane foam is highly durable and easy to clean

#58171 Rotating Arm, Inter-Lock Mount, Right or Left



Panacea® Foot – has a sleek, subtle design

- Injection-molded, closed-cell foam alleviates pressure
- Protects lower leg and foot from unwanted movement
- Heel cutout for better positioning, pressure relief
- Swings away with legrest for easy transfers
- Bolts secure to any footplate for impact protection
- Waterproof and easy to clean


#64647 Panacea Foot

#71600 Panacea Foot, pair

Panacea® Wheelchair Positioning Devices

Also Available

Panacea® Foam and Dynamic Air Mattresses

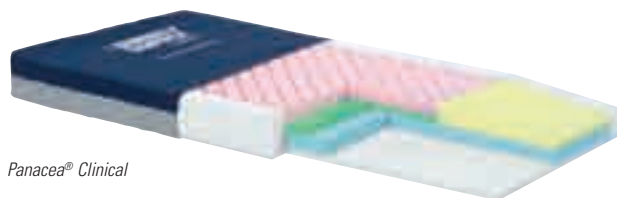

Panacea® Clinical



Panacea® Air Advance

Taking Cost Out of Healthcare®



directsupply.com

SS2706