

Direct Supply® Panacea® Air Max Bariatric Mattresses



Alternating pressure and true low air loss for residents up to 1,000 lbs.



#27505

Benefits

- Powerful support and performance for bariatric residents
- Three separate zones of support, including static head section, plus foot and mid zones
- Alternating pressure cycle times are adjustable from 10, 15, 20 or 25 minutes
- Powerful 1,000 lpm true low air loss wicks away heat and moisture
- Standard firm side perimeters protect against entrapment and minimize resident falls. Raised side perimeters available
- Each air cell is easily removable via quick disconnect for maintenance and replacement
- Air cells in head section do not alternate, for resident comfort
- Fowler mode provides extra support to resident when in the upright, seated position to prevent bottoming out
- Adjustable firmness levels for resident comfort and proper sizing
- Lockout mode prevents tampering and accidental adjustment

Specifications	
Weight Capacity	1,000 lbs. (42" – 54"W)
Therapy Style	Alternating Pressure
Firm Perimeters	Yes
Heel Slope	No
CPR Release	Yes
Safety Base*	Yes
Air Cell Height	8"
Mattress Depth	10"
Cover Type	4-way stretch
**Medicare HCPCS Code	E0277
Limited Warranty	Mattress – 2 years; Pump – 18 months

*Prevention against bottoming out during loss of power or air cell puncture

**This information is intended as a general guideline and is not billing or legal advice. Direct Supply does not guarantee coverage or reimbursement of any products. It is your responsibility to ensure the appropriateness of each claim submission, in accordance with payer requirements.

May be appropriate for up to Stage IV pressure injuries. Mattress wound Stage ratings are general usage guidelines. Resident-specific assessment could alter your particular usage of these mattresses.

SAVE THOUSANDS

by owning instead of renting!
See details at

[www.directsupply.com/blog/
renting-vs-owning-an-air-mattress-
in-long-term-care/](http://www.directsupply.com/blog/renting-vs-owning-an-air-mattress-in-long-term-care/)


SIZE		CONFIGURATION	
Width	Length	Flat	"Plus" Raised Perimeter
42"	80"	#27505	#L8338
	84"	#27506	#L8339
48"	80"	#18671	#18672
	84"	#2DM26	–
54"	80"	#32594	–



1-800-634-7328 / DirectSupply.com

Panacea® Air Max Bariatric Mattresses

Alternating pressure and true low air loss for residents up to 1,000 lbs.

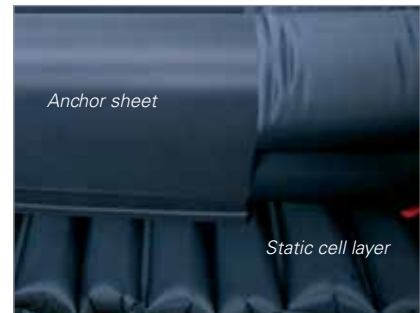

Mattress Body

- Twenty 8"H air cells provide superior envelopment and pressure redistribution for bariatric residents
- Coated air cells built to withstand bariatric use
- Each air cell is easily removable via quick-disconnect for maintenance
- Air cells are Velcro® attached to the top of side perimeters to minimize gaps and resident falls
- Firm side perimeters minimize the risk of falls and entrapment, help with resident transfers
- Optional 3"H raised side perimeters for additional fall prevention
- Anchor sheet holds air cells securely in the center to prevent them from separating or migrating from their intended location
- Side perimeters and additional bottom layer of 2"H cells are statically filled and remain inflated in the event of a power outage, for resident safety
- CPR quick-release valves allow for emergency deflation for resident safety

Cover

Low-shear, antimicrobial, fluid-resistant 70-denier cover

- 4-way stretch top cover wicks away heat and moisture, allows for maximum resident immersion
- Bottom cover: vinyl, fluid proof

Concealed zipper to prevent skin tears and protect foam from liquids

- Extra-large flap covers zipper to prevent contact with skin and prevent fluid from entering
- Zipper allows for easy cover removal and replacement if necessary

Control Unit

- Provides alternating pressure on a 5-, 10- or 15-minute cycle
- 1,275 lpm low air loss output for fast set-up and adjustments
- Low-pressure light and audible alarm indicate loss of pressure and power, increasing resident safety
- Fowler mode adjusts pressure in the mattress when resident is in the upright position to prevent bottoming out
- Low air loss only, alternating pressure with low air loss and max inflate modes
- Static mode aids in resident transfers and repositioning
- Adjustable firmness levels for resident comfort and proper sizing
- Lockout button prevents tampering and accidental adjustment
- Inflates in 50 seconds