
# Panacea® Original Bariatric Mattresses





This mattress may be appropriate for stage I or II pressure wounds.\*

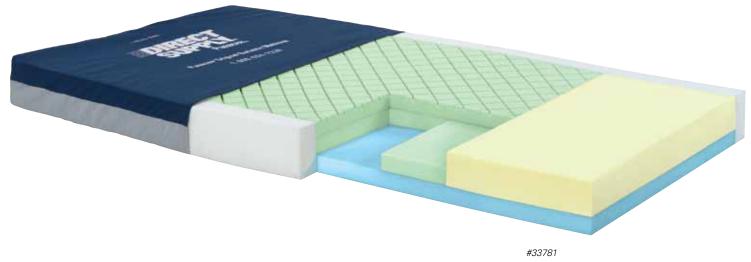

#### **Benefits**

- Die-cut, high-density top layer in the head and torso area provides comfortable pressure reduction
- Die-cuts create individual load-bearing cells to cradle residents, and allow air to flow and cool the sleep surface
- High-resiliency, 4"H super-soft foam heel area is softer than the mattress body to protect vulnerable heels
- Built-in firm perimeters make for easier transfers, and help reduce the risk of falls and entrapment
- Full-length base provides a firm, supportive foundation to prevent bottoming out
- 4-way stretch cover is low shear, antimicrobial, fluid resistant, breathable and allows for maximum resident immersion and envelopment; heavyduty, fluid-proof

Specifications		
Weight Capacity	800 lb.	
Standard Dimensions	42" or 48"W; 7"H; 80" or 84"L	
Heel Slope	No	
Raised Perimeters	No	
Cover Type	4-way stretch	
Meets 16 CFR 1633	Yes	
Limited Warranty	5-year, non-prorated limited warranty; 1-year for cover	


These pressure redistribution support surfaces are appropriate for use as part of an overall care plan to prevent and treat decubitus ulcers. Resident-specific assessment could alter your particular usage of these mattresses.


\*Mattress wound stage ratings are for guideline purposes only and are not claims for efficacy on individual residents. Mattresses are one part of an overall care plan to prevent and treat decubitus ulcers. A complete clinical assessment of each resident should be performed before selecting any therapeutic support surface.




## Direct Supply® Panacea® Original Bariatric Mattresses



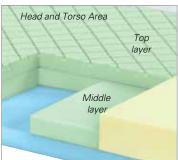

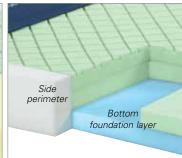

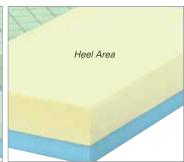

### Foam

#### Head and Torso Area Top Layer - 2"H

- Die-cut, soft foam top layer in the head and torso for effective pressure management
- Diamond-cut, die-cut top layer permits airflow to keep the mattress cool and dry as well as creates channels to reduce shear and friction

#### Head and Torso Middle Layer – 2"H

A firmer layer of foam to provide additional support for therapeutic top layer

#### Bottom Layer - 3"H

- A firmer layer of durable foam, providing support and longevity for the entire mattress
- · Prevents resident from bottoming out

#### Heel Area - 4"H Top Layer

- Ultra soft, high-resiliency foam to cradle the vulnerable heel area Firm Side Perimeters
- Built-in firm foam perimeters make for easier transfers and help prevent falls

#### Cover

Low-shear, antimicrobial, fluid-resistant

- Four-way stretch top cover wicks away heat and moisture, allows for maximum resident immersion
- Bottom cover: vinyl, fluidproof

Concealed zipper to prevent skin tears and protect foam from liquids

- Extra-large flap covers zipper to prevent contact with skin and prevent fluid from seeping in
- Zipper allows for easy cover removal and replacement if necessary

Width	Length	Flat
42"	80"	#33781
48"		#33782
42"	84"	#33789
48"		#33790

