

# MATTRESS LATERAL ROTATION



## COMBINES 3 THERAPIES

- LATERAL ROTATION
- ALTERNATING PRESSURE
- LOW AIR LOSS  
TO EFFECTIVELY TREAT PRESSURE ULCERS.



### LATERAL ROTATION

Rotates the patient from side to side to provide proper pressure redistribution.



### ALTERNATING PRESSURE

Redistributes pressure throughout the length of mattress.



### LOW AIR LOSS

Treats bed sores by creating a dry cool environment.





Customized  
Lateral Rotation



Side bolsters for optimum patient safety reduces the risk of falling.

Mid-way transfer point allows for easy patient egress and ingress

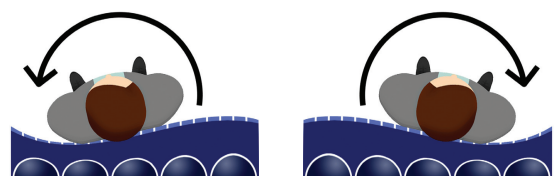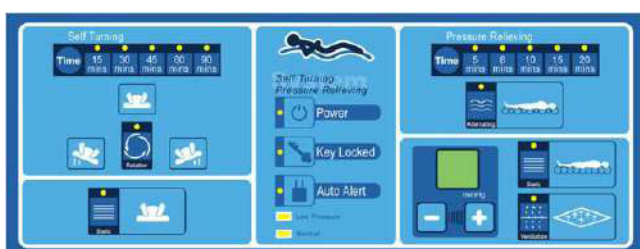

# MATTRESS LATERAL ROTATION



## PUMP



Alternating Pressure and Static Function



Comfort Level Adjustable



Low Air Loss with on/off function



5 Level Cycle Time Adjustable



- 5 Level Rotation Time Adjustable
- Rotation On/Off Function
- Panel Lock Function
- Power Failure Alert
- Low Pressure Alert
- Audible Alarm Mute Function
- Quick-change Air Filter

### PRESSURE RANGE

Digital readout adjusts pressure range from 30 - 60mmhg

### ALTERNATING CYCLE TIME

5, 8, 10, 15, 20 minutes selectable

### ROTATION CYCLE TIME

15, 30, 45, 60, 90 minutes selectable

### AIR OUTPUT

10 LPM

### PUMP DIMENSIONS

11" x 12" x 5.9"

## MATTRESS



Alternating Pressure



Low Air Loss



Lateral Rotation



Pillow Function - Top 3 cells remain static for ultimate patient comfort



Adjustable Straps secures mattress to bed frame



### LOW AIR LOSS

8 Ventilated Cells

### CELL NUMBER

20 Easy Replaceable Air Cells

### CELL MATERIAL

TPU

### COVER MATERIAL

2 Way Stretch PU Fluid Resistant Cover

### BASE MATERIAL

Woven with coated non-slip material

### CPR

Pull strap

### DIMENSIONS

80" L x 36" W x 8" H with 2" Side Bolsters

### WEIGHT CAPACITY

440 Lbs.



Safe yet effective Angle of Rotation: 30 degrees



MODEL # **LRM36**