

# Direct Supply® Panacea® Air Max Bariatric Mattresses



Alternating pressure and true low air loss for residents up to 1,000 lbs.



## Benefits

- Powerful support and performance for bariatric residents
- Three separate zones of support, including static head section, plus foot and mid zones
- Alternating pressure cycle times are adjustable from 10, 15, 20 or 25 minutes
- Powerful 1,000 lpm true low air loss wicks away heat and moisture
- Standard firm side perimeters protect against entrapment and minimize resident falls. Raised side perimeters available
- Each air cell is easily removable via quick disconnect for maintenance and replacement
- Air cells in head section do not alternate, for resident comfort
- Fowler mode provides extra support to resident when in the upright, seated position to prevent bottoming out
- Adjustable firmness levels for resident comfort and proper sizing
- Lockout mode prevents tampering and accidental adjustment

Specifications	
<b>Weight Capacity</b>	1,000 lbs. (42" – 54"W)
<b>Therapy Style</b>	Alternating Pressure
<b>Firm Perimeters</b>	Yes
<b>Heel Slope</b>	No
<b>CPR Release</b>	Yes
<b>Safety Base*</b>	Yes
<b>Air Cell Height</b>	8"
<b>Mattress Depth</b>	10"
<b>Cover Type</b>	4-way stretch
<b>Limited Warranty</b>	Mattress – 2 years; Pump – 18 months

\*Prevention against bottoming out during loss of power or air cell puncture

*May be appropriate for up to Stage IV pressure injuries. Mattress wound Stage ratings are general usage guidelines. Resident-specific assessment could alter your particular usage of these mattresses.*

## SAVE THOUSANDS

by owning instead of renting!  
See details at

[www.directsupply.com/blog/  
renting-vs-owning-an-air-mattress-  
in-long-term-care/](http://www.directsupply.com/blog/renting-vs-owning-an-air-mattress-in-long-term-care/)


SIZE		CONFIGURATION	
Width	Length	Flat	"Plus" Raised Perimeter
42"	80"	#27505	#L8338
	84"	#27506	#L8339
48"	80"	#18671	#18672
	84"	#2DM26	–
54"	80"	#32594	–



1-800-634-7328 / [DirectSupply.com](http://DirectSupply.com)

# Panacea® Air Max Bariatric Mattresses

Alternating pressure and true low air loss for residents up to 1,000 lbs.

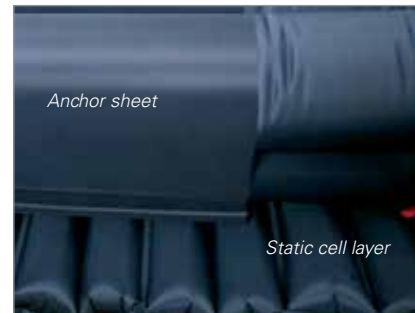

## Mattress Body

- Twenty 8"H air cells provide superior envelopment and pressure redistribution for bariatric residents
- Coated air cells built to withstand bariatric use
- Each air cell is easily removable via quick-disconnect for maintenance
- Air cells are Velcro® attached to the top of side perimeters to minimize gaps and resident falls
- Firm side perimeters minimize the risk of falls and entrapment, help with resident transfers
- Optional 3"H raised side perimeters for additional fall prevention
- Anchor sheet holds air cells securely in the center to prevent them from separating or migrating from their intended location
- Side perimeters and additional bottom layer of 2"H cells are statically filled and remain inflated in the event of a power outage, for resident safety
- CPR quick-release valves allow for emergency deflation for resident safety

## Cover

Low-shear, antimicrobial, fluid-resistant 70-denier cover

- 4-way stretch top cover wicks away heat and moisture, allows for maximum resident immersion
- Bottom cover: vinyl, fluid proof

Concealed zipper to prevent skin tears and protect foam from liquids

- Extra-large flap covers zipper to prevent contact with skin and prevent fluid from entering
- Zipper allows for easy cover removal and replacement if necessary

## Control Unit

- Provides alternating pressure on a 5-, 10- or 15-minute cycle
- 1,275 lpm low air loss output for fast set-up and adjustments
- Low-pressure light and audible alarm indicate loss of pressure and power, increasing resident safety
- Fowler mode adjusts pressure in the mattress when resident is in the upright position to prevent bottoming out
- Low air loss only, alternating pressure with low air loss and max inflate modes
- Static mode aids in resident transfers and repositioning
- Adjustable firmness levels for resident comfort and proper sizing
- Lockout button prevents tampering and accidental adjustment
- Inflates in 50 seconds